

on/go™ COVID-19 Antigen Self-Test

Rapid Diagnostic Test for Detection of SARS-CoV-2 Antigen


Instructions for Running the Test

IMPORTANT: Do not open kit components until instructions to do so. Follow the in-app self-paced, step-by-step instructions or paper instructions printed on the QRI as below.

- 1** Wash your hands thoroughly for at least 20 seconds before the test.



- 2** Unpack the test components from the tray.

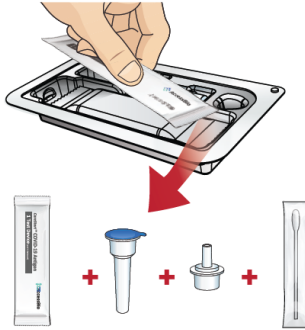


- 3** Remove the test cassette and place it on a flat, clean surface.




- 4** Locate the extraction vial and gently peel off the aluminum foil seal, being sure to keep the vial upright and place it in the packaging tray.



- 5** Locate a nasal swab and remove from the pouch. Be careful not to touch the swab tip.

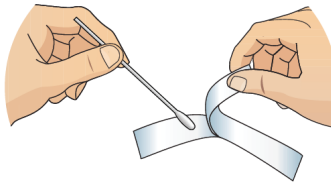


- 6** Gently insert the swab no more than 3/4 inch into the **LEFT** nostril. Then, slowly rotate the swab at least **5 times** in a circular path for a total of **15 seconds**. If you have questions, see the CDC Guidelines.

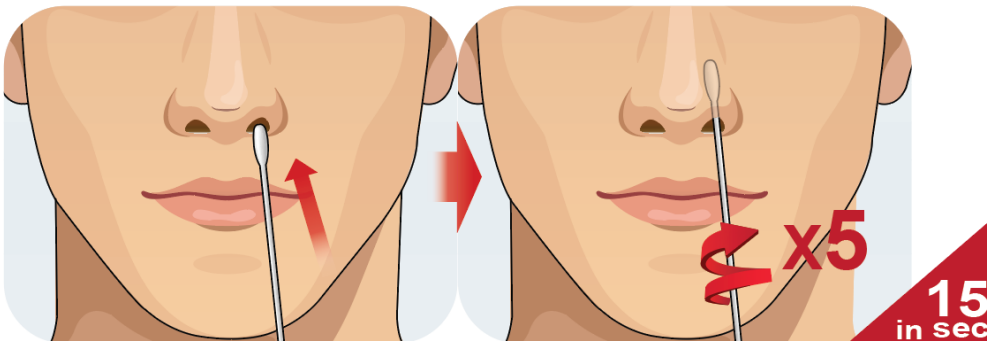


on/go™ COVID-19 Antigen Self-Test

Rapid Diagnostic Test for Detection of SARS-CoV-2 Antigen


- 7** Gently remove the swab from the LEFT nostril and place directly into the **RIGHT** nostril, repeating the process of rotating at least **5 times** in a circular path for a total of **15 seconds**. Remove the swab from the RIGHT nostril.

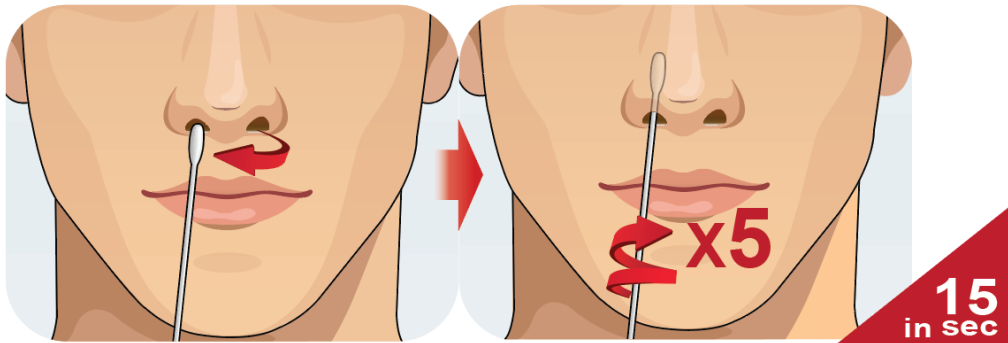


- 8** Place the swab into the extraction vial. Rotate the swab vigorously at least 5 times.

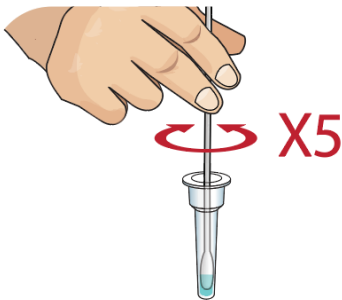

- 9** Remove the swab by rotating against the extraction vial while squeezing the sides of the vial to release the liquid from the swab. Discard the swab in trash.

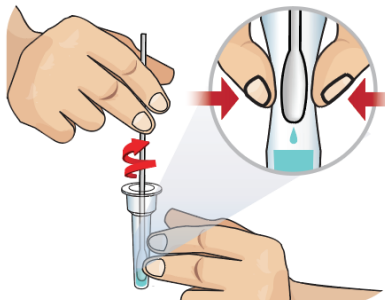


- 10** Close the vial by pushing the cap firmly onto the vial.

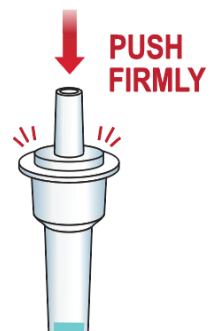

- 11** With your finger, mix thoroughly by flicking the bottom of the vial.



- 12** Invert the extraction vial and hold the sample vertically above the sample well. Squeeze the vial gently. Allow **THREE (3)** drops of sample to fall into the sample well.

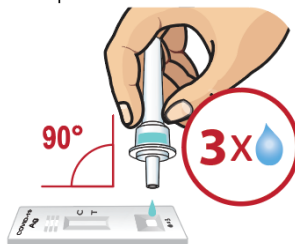

- 13 Start a timer.**


Read the result at **10 minutes**.
The test result should not be read after 15 minutes.

**IMPORTANT**

Do not move or lift the test cassette during this time.
Do not exit the mobile app during this process.

Disposal

Dispose of all used test kit components and swab samples in household trash.

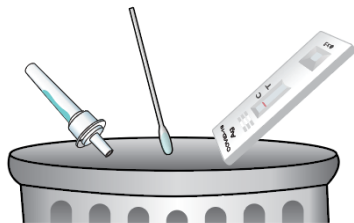

If used for serial testing and the test result is negative, a second test should be obtained two or three days with at least 24 hours and no more than 48 hours between tests.


Interpretation of Results

The test results will be interpreted by visual reading following the in-app interpretation instructions or provided Quick Reference Instructions.

NOTE: The test results should be read by visual and interpreted at 10 minutes after the sample application and the reading and interpretation of the results should not exceed 15 minutes as it may yield inaccurate results.

COVID-19 Detected (Positive):

One purple-colored line next to “C” and one blue-colored line next to “T” indicates COVID-19 positive result.

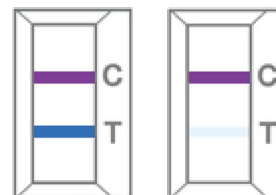


NOTE: The color intensity of the blue-colored test line will vary depending on the amount of SARS-CoV-2 nucleocapsid protein antigen present in the sample. Any faint blue-colored line in the test line should be considered as positive.

A positive test result for COVID-19 indicates that antigens from SARS-CoV-2 were detected, and the patient is very likely to be infected with the virus and presumed to be contagious. Test results should always be considered in the context of clinical observations and epidemiological data (such as local prevalence rates and current outbreak/epicenter locations) in making a final diagnosis and patient management decisions. Patient management should follow current CDC guidelines. Additional confirmatory testing with a molecular

test for positive results may also be necessary, if there is a low likelihood of COVID-19, such as in individuals without known exposures to COVID-19 or residing in communities with low prevalence of infection.

COVID-19 Not Detected (Negative):

One purple-colored line only next to “C” indicates a negative result.

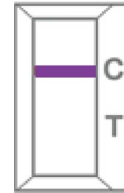


Negative results do not rule out SARS-CoV-2 infection. Individuals without symptoms that test negative should be tested again with at least 24 hours and no more than 48 hours between tests. All negative results are considered presumptive, and confirmation with a molecular assay, if necessary for patient management, may be performed. The amount of antigen in a sample may decrease as the duration of illness increases. Negative results should be considered in the context of an individual's recent exposures, history, and the presence of clinical signs and symptoms consistent with COVID19 and confirmed with a molecular assay, if necessary, for patient management.

Invalid:

Invalid barcode or absence of a purple-colored line next to “C”. Re-test with a COVID-19 test may be needed.

An invalid test result indicates that your test has experienced an error and unable to interpret the result of the test. You will need to retest with a new test or consult a healthcare professional. If you still have symptoms, you should self-isolate at home and avoid contact with others prior to the retest

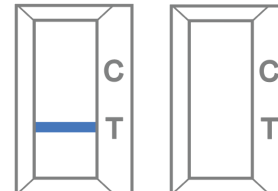

For questions, or to report a problem, please call Technical Support at +1-888-965-0301 or visit letsongo.com/Support.

Limitations

1. This test will indicate the presence of SARS-CoV-2 nucleocapsid protein antigen in the specimen from both viable and non-viable SARS-CoV-2 virus. Test performance depends on the amount of virus (antigen) in the sample and may or may not correlate with viral culture results performed on the same sample.